

EPILEPSY & The Facts

Vagal Nerve Stimulator (VNS)


What is Vagus Nerve Stimulation?

Vagus Nerve Stimulation (VNS) is the mild electrical stimulation of a nerve called the vagus nerve that runs through the neck and carries information to the brain. A pacemaker-like generator implanted just under the skin in the patient's left chest area delivers electrical impulses to this nerve at regular intervals around the clock, 24 hours a day, every day.

How can Vagus Nerve Stimulation control seizures?

It is thought that nerve fibers from the vagus nerve carry information to the brain stem and thence to areas of the brain known to be involved in seizures. Epileptic seizures are caused by normally well-organized brain waves becoming completely disorganized. It has been suggested that stimulation of the vagus nerve may be able to interrupt this abnormal brain activity, reducing or possibly preventing seizures.

What does the Vagus Nerve Stimulation procedure involve?

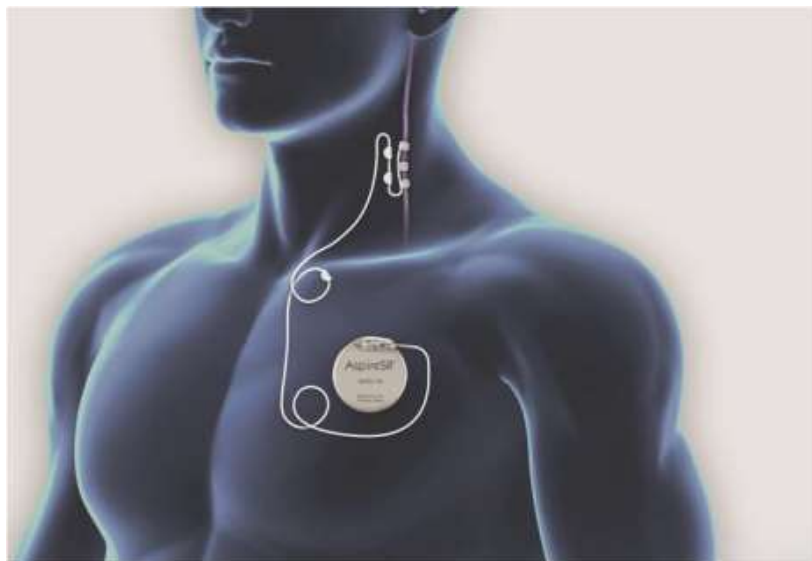

The VNS Therapy system is comprised of two parts: the generator and the lead (wire that connects the generator to the vagus nerve). Typically surgery requires 60 – 90 minutes, and is performed under general anesthesia.

- Two incisions are made, each about 5-6 cm long, one on the skin in the left chest area, and the other on the left side of the neck.
- A small pulse generator, about the size of a stopwatch, is surgically implanted in the patient's upper left chest area.
- A wire is tunneled under the skin and attached to the vagus nerve in the neck.
- The electrical connection between the vagus nerve and the pulse generator is checked before the incision is closed.

Either at time of implantation or after a short healing period, the battery-powered pulse generator is programmed by the doctor to begin sending a particular “dosage” of stimulation to the vagus nerve, typically 30 seconds of stimulation, at five-minute intervals, with a strength of 0.25 milliamperes (mA) of electrical current. Initially, more frequent follow-up visits are needed to gradually increase the output of higher levels of stimulation in order to achieve better levels of seizure control. Once the ideal “recipe” for stimulation, combining device output and on/off sequencing, has been determined for the patient, follow-up with the neurologist returns to usual frequency.

Can the patient activate their Vagus Nerve Stimulator?

Some people with epilepsy experience an aura (which is really a small seizure) warning them that a bigger seizure is about to occur. In these cases, passing a hand-held magnet over the pulse generator can activate it. This sends an extra burst of stimulation to the vagus nerve, which may stop the seizure or make it less severe. A caregiver or companion can also use the magnet on behalf of the patient.

Unfortunately, the timing or circumstances of the seizure result in a situation where the magnet cannot always be used, or is not always available. People who experience night-time (nocturnal) seizures are unlikely to awaken prior to the seizure to be able to use the magnet. For some people, the seizure onset is so sudden that there is insufficient time to apply the magnet. Still others, due to physical limitations or the effect of the seizure itself, simply cannot physically apply the magnet. Note that a caregiver can apply the magnet, and even using the magnet partway into a seizure may have benefit.

Nevertheless, for these reasons, the newest version of the VNS Therapy device uses a surrogate for seizure activity, and automatically triggers an extra dose of stimulation. Most patients experience a sudden heart rate increase at the time of at least one of their seizure types. This heart rate rise (called an ‘ictal tachycardia’ [IK-tal ta-kee-KAR-dee-ah]) can be detected by the implanted system. It looks for a sudden jump in heart rate (not the *gradual* heart rate increase experienced in exercise), and then automatically provides an extra dose of stimulation. Thus, it can be considered a ‘magnet-less’ magnet mode. The settings for this function, known as ‘autostimulation’, can be tailored to the patient’s individual heart rate response and needs. The autostimulation mode is used in conjunction with the normal operation and also the magnet function.

What are the side effects of Vagus Nerve Stimulation?

As with any surgery, there is a risk of infection. When the vagus nerve is actually being stimulated, many people experience a tingling in the neck or hoarseness in their voice. Some have reported throat, ear or tooth pain, coughing and shortness of breath. Having the intensity or duration of the stimulation adjusted by the doctor may reduce or stop these sensations. Importantly, most people report a reduction in side effects over time, such that by three years post-implantation very few people report any issues at all.

Because VNS Therapy has been used for over 20 years now, the side effect profile has been well established in multiple scientific publications. VNS Therapy is safe for long-term use, and pregnancies have gone to term in women receiving VNS Therapy.



The Critical Care Services Ontario (CCSO) publication in March 2016 of "Provincial Guidelines for the Management of Medically-Refractory Epilepsy in Adults and Children who are not Candidates for Epilepsy Surgery" (available on-line at

http://braininstitute.ca/docs/provincial_guidelines_for_management_of_mre_in_adults_and_children_who_are_not_candidates_for_surgery.pdf) has identified the need for specific funding for epilepsy surgery that is apart from other neurosurgeries. This program has been implemented, and four centres have been identified as Regional Epilepsy Surgery Centre of Excellence (RESC) as shown in the grid below:

	Pediatric	Adult
Toronto	The Hospital for Sick Children (SickKids)	Toronto Western Hospital
London	London Health Sciences Centre (Victoria Hospital)	London Health Sciences Centre (University Hospital)

Fully funded VNS Therapy implantation surgery is performed at these hospitals. There is no cost to the individual as it is covered under OHIP. In addition to these centres, routine follow-up of VNS Therapy devices (clinic) is available at these centres:

	Pediatric	Adult
Ottawa	Children's Hospital of Eastern Ontario (CHEO)	Ottawa General Hospital
Kingston	Hotel-Dieu Hospital	Kingston General Hospital
Hamilton	McMaster University Hospital	Hamilton General Hospital

Other clinics/physicians may also provide follow-up service beyond these programmes.

